

*Otto Bock®*


## Kimba Cross



A harmonious blend of function and design. Oval-crossection, aluminium tubing and a telescoping frame make the Kimba Cross highly functional: for example, it adapts to fit various-sized seating units. Furthermore, the Kimba Cross incorporates the Kimba seating system, permitting a wide range of positioning possibilities and individual support. The adjustable wheel base permits high stability.

Thanks to a low centre of gravity, a spring system for the seat, large pneumatic tires, and a drum brake, Kimba Cross offers a high level of security in demanding everyday situations, for easy handling even in rugged terrain. Whether indoors or out: the swivelling front wheel or the cycle trailer package extend the radius of action. – Kimba Cross: attractive and functional.

**The All-Rounder for Everyday and Leisure.**

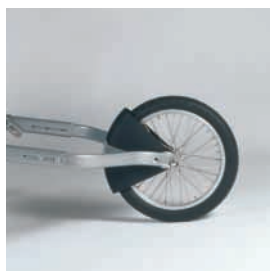

## Kimba Cross in Detail

### Outdoor Mobility Base

- Telescoping frame fits various-sized seating units
- Spring system for the seat
- 16" front wheel, pneumatic tire, rigid, quick-lock attachment, incl. splash guard
- 18" rear wheels, quick-release axle
- Drum brake with lock option
- Height-adjustable pushbar
- Can be folded and disassembled for easy stowage
- Attractive design thanks to oval tubular frame

### Kimba Multifunctional Seating Unit

- Optimal adjustability to body proportions, individual positioning thanks to support pad system and fixation possibilities, numerous possibilities for adjustment / variation of seat position

### Options

- Swivelling wheel with quick-lock attachment for confined spaces
- Cycle trailer package (draw-bar with quick-lock, cycle coupling, pennant, light)
- Anti-tipper, clothing protectors
- Spoke protectors, storage bag

### Technical Data

	Size 1: 1-6 years	Size 2: 4-10 years		
Dimensions in cm/kg	Outdoor mobility base		Seat Size 1	Seat Size 2
Pushbar height	90-115			
Width	66			
Length	140-170			
Height	75-120			
Seat depth			19-31	24-40
Seat width			20-30	26-35
Lower leg length			19-31	20-37
Backrest height/angle			46-66/0° to +90°	56-70/0° to +90°
Wheel diameter front/back	40/45 (swivelling wheel 30/45)			
Seat tilt	-5° to +30°	-		
Load capacity	50		35	35
Weight	15.8 (incl. splash guard)		7	8.8
Folding size (LxWxH)	94x60x32		72x40x22	82x44x27