




Increasingly therapists and families are realising the benefits to be gained from riding tricycles. Not only are they seeing the fun triking gives the rider they are witnessing the therapeutic benefits it provides through muscle stimulation and cardiovascular exercise.

Since the first tricycle was conceived and produced in 1970, we have sought to design tricycles which are appealing to children and adults of all ages. Whether for use at home, at school or for cycle expeditions, our trikes are designed for everyone to enjoy. The secure feeling provided by the broad wheelbase and well balanced handling coupled with the extensive range of positional accessories ensures that they can accommodate most physically challenged riders.

Available in five sizes; Imp, Terrier, TMX, Tracker 20" and Tracker 24", these trikes share the same distinctive sporty appearance. Comfort and support is provided by the extra wide padded saddle, an essential component therapists recommend for a stable seating position. This is complemented by the adjustable High-Low bend handlebars and seat adjustments which combine to ensure the rider is in the optimum position to pedal and control the tricycle in a confident manner.


## Accessories

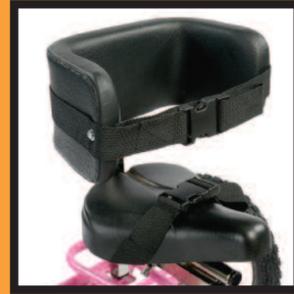

Adjustable back/trunk support - Imp



Adjustable back/trunk support - Terrier, TMX, Tracker

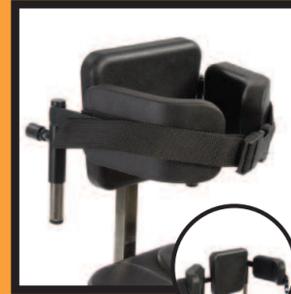


In-line trunk support




Comfort saddle




Footsandals C/W pulley system

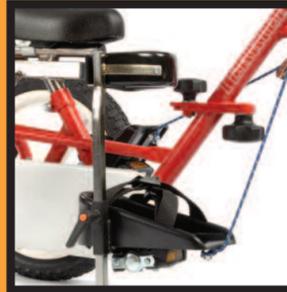


Leg calipers

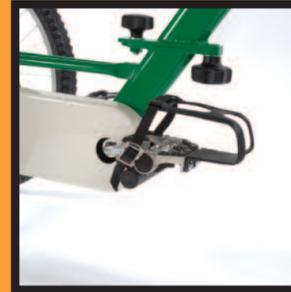


Self balance pedals C/W toe clips

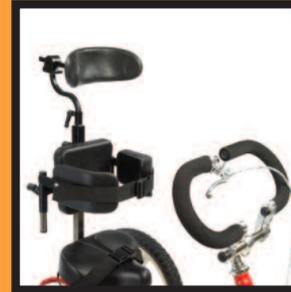

Adjustable head support

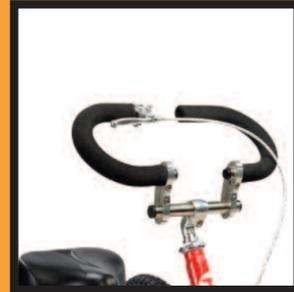

Cumfigrip handlebar

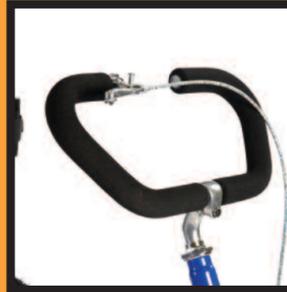

Loop handlebar



Push/pull pole



Rear steering system



# Tricycles



32 Welbeck Road, Darnley Industrial Estate, Glasgow, G53 7SD, Scotland  
Tel: 44 (0) 141 876 9177 Fax: 44 (0) 141 876 9039


Visit us at: [www.theraplay.co.uk](http://www.theraplay.co.uk)  
PO Box 1364, Cumberland, MD 21501-1364, USA  
Tel: 301 759 3525 or 800 306 6777 Fax: 301 759 3525

Visit us at: [www.triaid.com](http://www.triaid.com)

**A whole new world of possibilities**

# Imp


Designed with younger children in mind, it is suitable for children from 2½ years. It has a low gear ratio for easy pedaling, 12½" wheels, parking brake and smaller handlebars for those smaller fingers and hands.

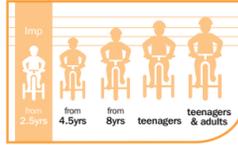

Model shown with adjustable trunk support and foot sandals

# Terrier

A rugged trike for children from 4½ years, this model has 16" wheels and a caliper brake with a parking brake mechanism. The terrier has a frame folding mechanism option which allows the frame to be folded for transportation and storage.

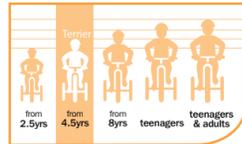

Model shown with adjustable trunk support and foot sandals



## Put play into *Therapy*

### Appearance

Looking good and feeling comfortable is important. Modern styling, the use of quality cycle components and a choice of ten frame colours mean that we have something for everyone.



### Frame Geometry

Designed to make pedaling easy, the frame geometry enables the rider to have a smooth even pedaling motion while maintaining a comfortable position.

### Gearing Options

All models are fitted with a low ratio 1:1 fixed gear (1 revolution of the pedals :1 revolution of the wheels). This enables riders to get going easily, allows the tricycle to be reversed and helps to slow down speeders. There are options to increase or decrease the gear ratio depending upon the riders ability. A derailleur gear option is available on the TMX and Tracker models.



### Dual Axle Drive System

A unique dual axle drive system has been developed to reduce resistance when pedaling. This feature ensures that our trikes are the easiest to pedal. If used in conjunction with a low gearing option, riders with very low muscle tone can enjoy the freedom of cycling independently.



### Stability

All models are designed with a low centre of gravity, a broad wheelbase and a wide rear axle for stability. A steering stop limiter is fitted to control the turning circle of the tricycle.

### Optional Accessories

A range of interchangeable accessories enable customization to suit a riders positional and support requirements.

### Colours

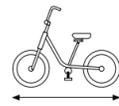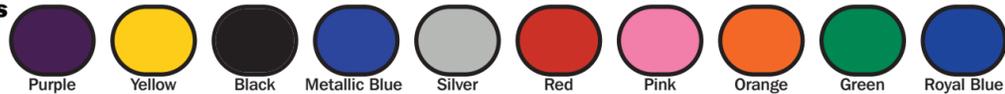

	Overall Length	Overall Width	To suit inside leg		Seat to handlebars		Max user weight
			(MIN)	(MAX)	(MIN)	(MAX)	
Imp	39½" (1000mm)	25½" (650mm)	12½" (320mm)	15½" (395mm)	11" (280mm)	14½" (370mm)	70lbs (32kg)
Terrier	45" (1140mm)	29" (740mm)	16½" (420mm)	21½" (545mm)	14" (360mm)	17" (430mm)	100lbs (45kg)
TMX	59" (1500mm)	29" (740mm)	20" (510mm)	25" (640mm)	19" (480mm)	22" (560mm)	130lbs (64kg)
Tracker 20"	65" (1650mm)	32" (810mm)	24" (610mm)	29" (740mm)	21" (535mm)	24" (610mm)	170lbs (77kg)
Tracker 24"	69" (1753mm)	32" (810mm)	25" (640mm)	30" (762mm)	22" (560mm)	25" (640mm)	170lbs (77kg)

# TMX


A robust tricycle designed to accommodate the larger and heavier child. The TMX has 20" wheels, a caliper brake with a parking brake mechanism and has a frame folding mechanism option which allows the frame to be folded for transportation and storage.



Model shown with in-line trunk support and foot sandals

# Tracker

The Tracker model is available in two sizes, the Tracker 20" is designed for teenagers and the Tracker 24" for teenagers and adults. Both models have alloy wheels, a caliper brake with a parking brake mechanism and are fitted with a frame folding mechanism which allows easy transportation and storage.

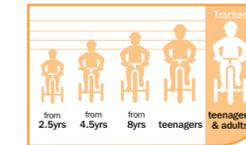

Model shown with folding option, tractor saddle and self balance pedals